

Using Illumina Sequencing Technology to Explore The Bacterial Genome

Björn Espedido, PhD

*Senior Field Application Scientist,
Illumina*

Zhiliang Chen, PhD

*Senior Bioinformatics Sales Specialist,
Illumina*

Agenda

SECTIONS



BaseSpace[®]
SEQUENCE HUB

1

Understanding the Illumina MiSeq system and library preparation

2

Introduction to BaseSpace Sequence Hub

3

Understanding bacterial whole genome alignment and variant calling


MiSeq and Library Preparation






MiSeq System Overview

MiSeq System Features

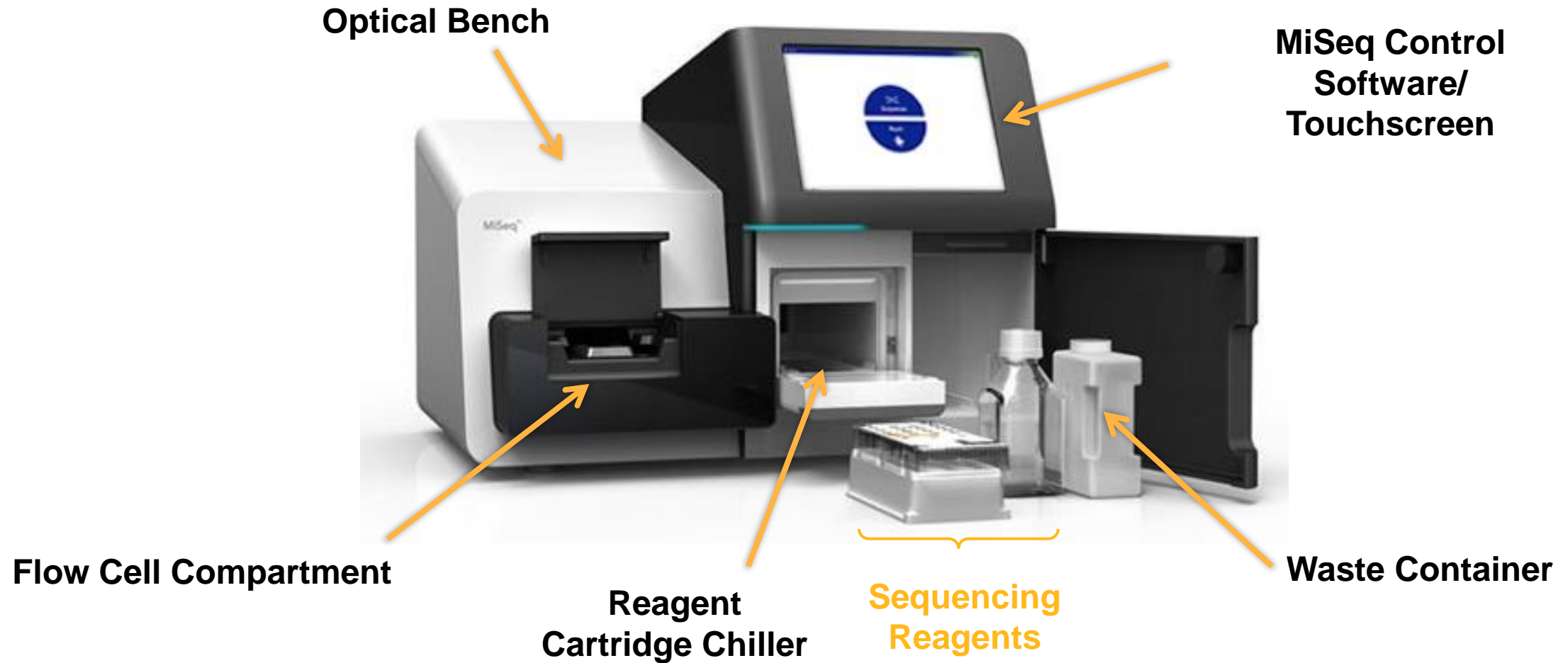

MiSeq Software


MiSeq Control Software (MCS)

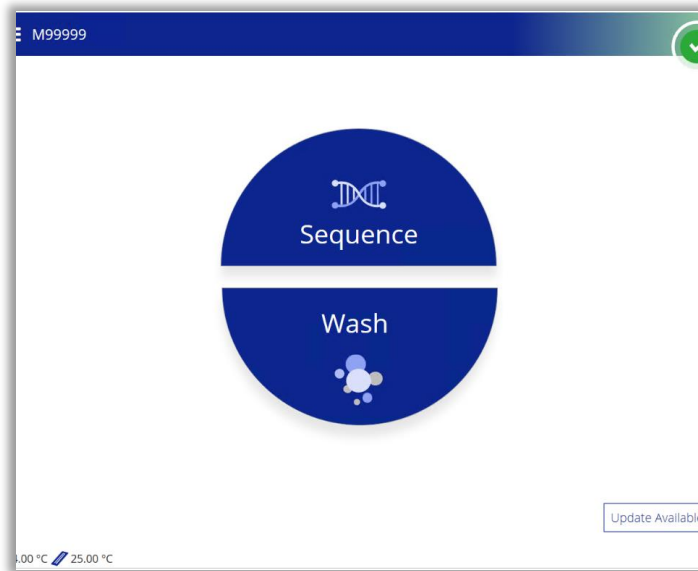

Main functions

- System operation
- Performing sequencing runs
- Performing primary analysis

Local Run Manager (LRM)

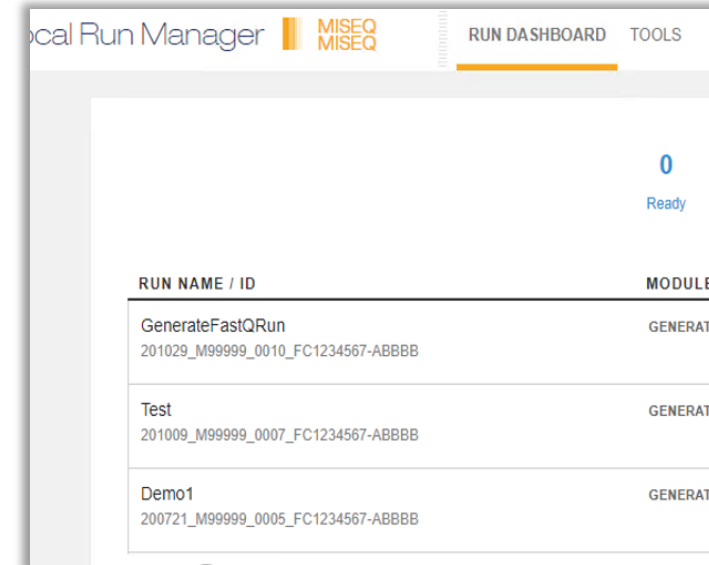

Main functions

- Creating a run
- Run monitoring
- Performing Secondary analysis


MiSeq Sequencing Reagent Kit

Four components, all single-use



Box	Components	Storage
Box 1/2	① Reagent Cartridge	Freezer, -25°C to -15°C
	② HT1	
Box 2/2	③ Flow Cell	Refrigerator, 2°C to 8°C
	④ PR2 Buffer	

Flow Cell Types

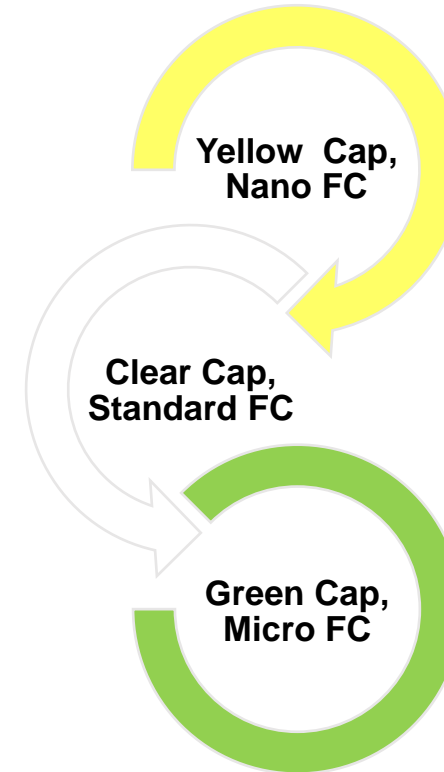

Use the reagent cartridge associated with your flow cell type

MiSeq Reagent Kit Variations

Standard Flow Cells



MISEQ REAGENT KIT V2

READ LENGTH	OUTPUT
1 × 36 bp	540-610 Mb
2 × 25 bp	750-850 Mb
2 × 150 bp	4.5-5.1 Gb
2 × 250 bp	7.5-8.5 Gb

MISEQ REAGENT KIT V3

READ LENGTH	OUTPUT
2 × 75 bp	3.3-3.8 Gb
2 × 300 bp	13.2-15 Gb

Reads Passing Filter	V2	V3
Single reads	12-15 M	22-25 M
Paired- End Reads	24-30 M	44-50 M

Low Output Applications


Flow cell	# of Reads	Read length	2 x 75 Output	2 x 150 Output	2 x 250 Output
Micro FC	4 M	Up to 2 x 150	600 Mb	1.2 Gb	----
Nano FC	1 M	Up to 2 x 250	150 Mb	300 Mb	500 Mb

Library Prep

Illumina Sequencing Workflow


1

Library Preparation

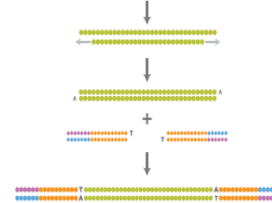


2

Cluster Generation

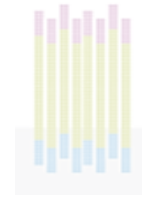

3

Sequencing

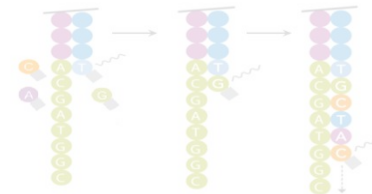


4

Data Analysis

```
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
```


Planning Your Prep and Run

- Before proceeding to library prep and creating your run, you will need to plan your run by defining the following parameters of your sequencing run:

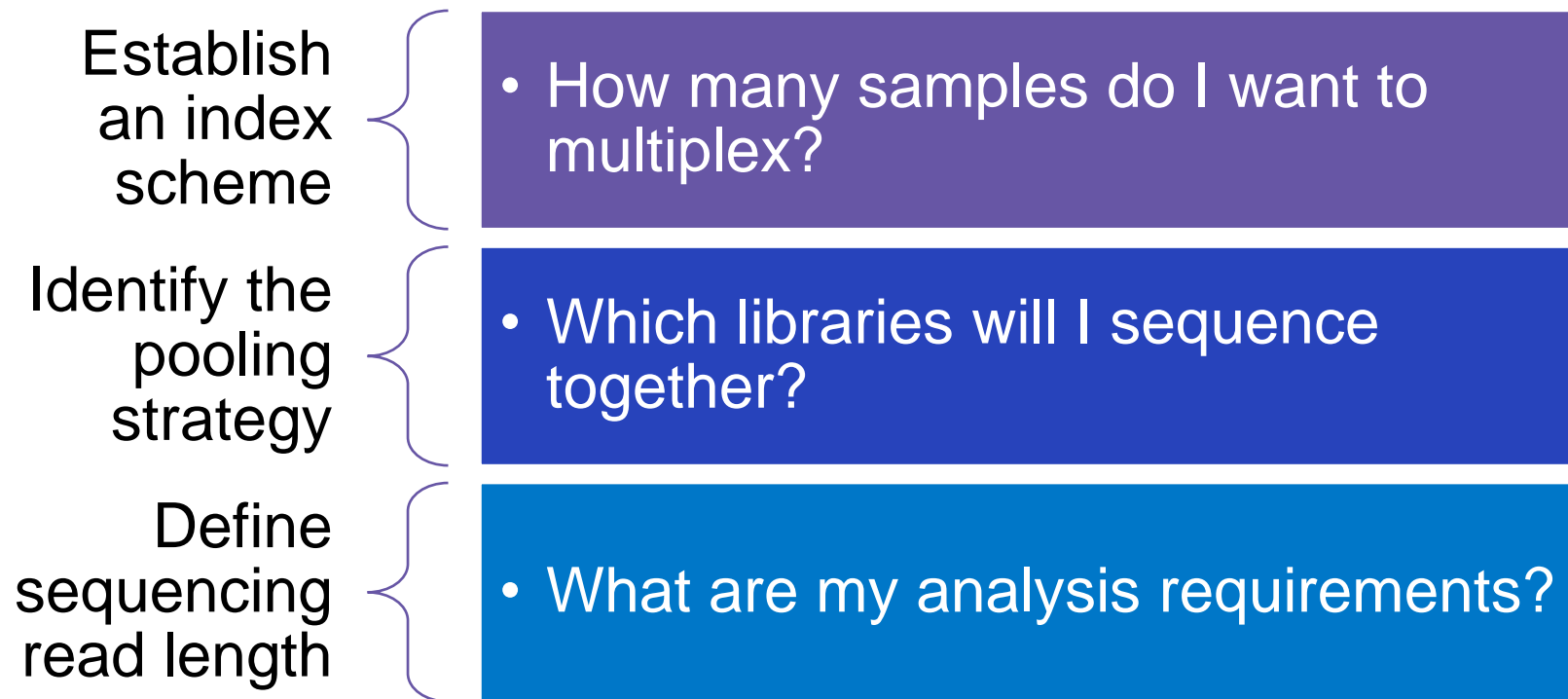

Library Prep is Critical for Successful Sequencing

The aim of library prep is to obtain nucleic acid fragments with adapters attached on both ends

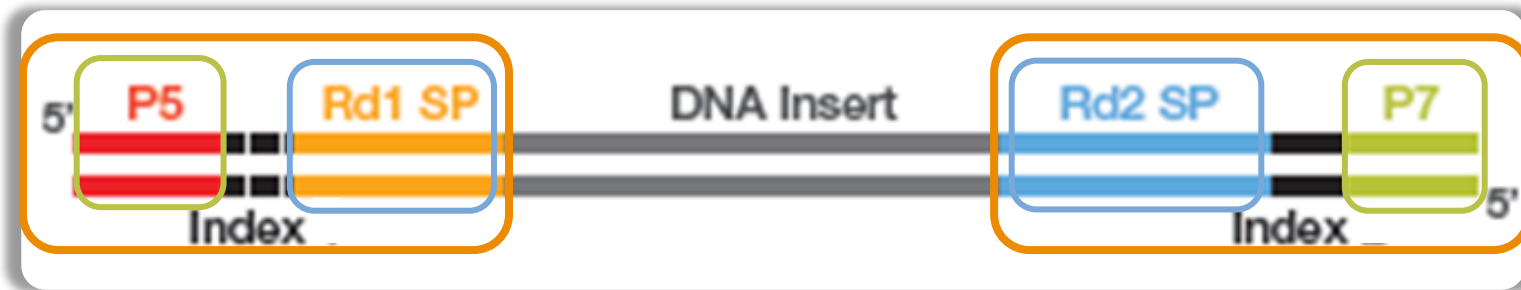

Dual Index Library shown

P5 and P7 sequences are complimentary to oligos bound to flow cell surface and *required* for any library




Indexes are used to tag individual samples to allow pooling of multiple libraries

Rd1 and Rd2 sequencing primers regions are used to initiate sequencing

Library Prep

Illumina DNA Prep, Tagmentation and DNA Indexes

Illumina Bead-Linked Transposome Chemistry

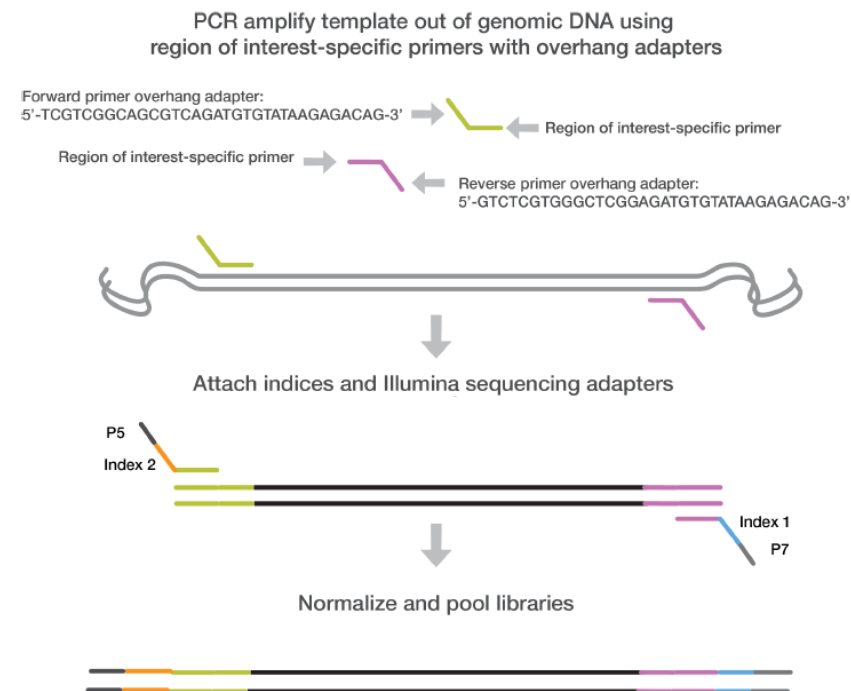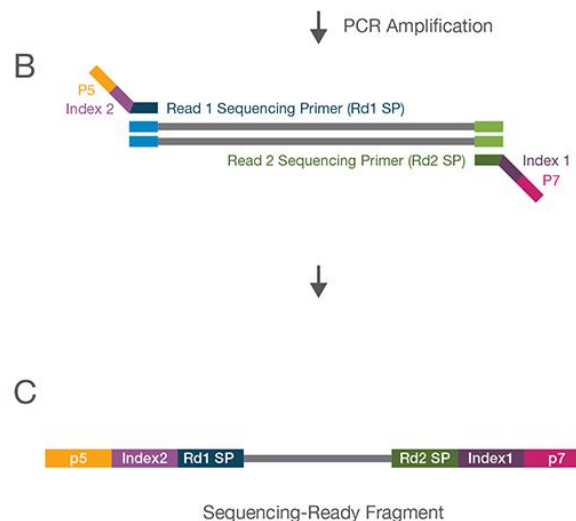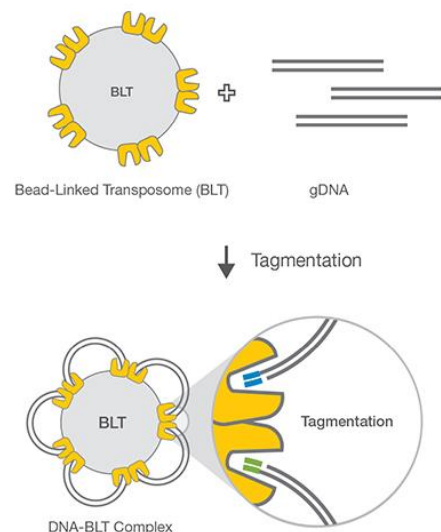

For large amplicons & whole genome sequencing

Amplicon sequencing using custom locus-specific primers with Illumina DNA indexes

Library Quality and Quantity Methods

- How are you planning to assess library quality and quantity?
 - Accurate quantification of quality libraries is critical for optimal cluster density and data quality

Recommended Quantity Methods

qPCR

Qubit

PicoGreen



Recommended Quality Methods

BioAnalyzer

Fragment Analyzer

Not Sure Which Kit to Use?

Use the Sequencing Sample Preparation Kit Selector

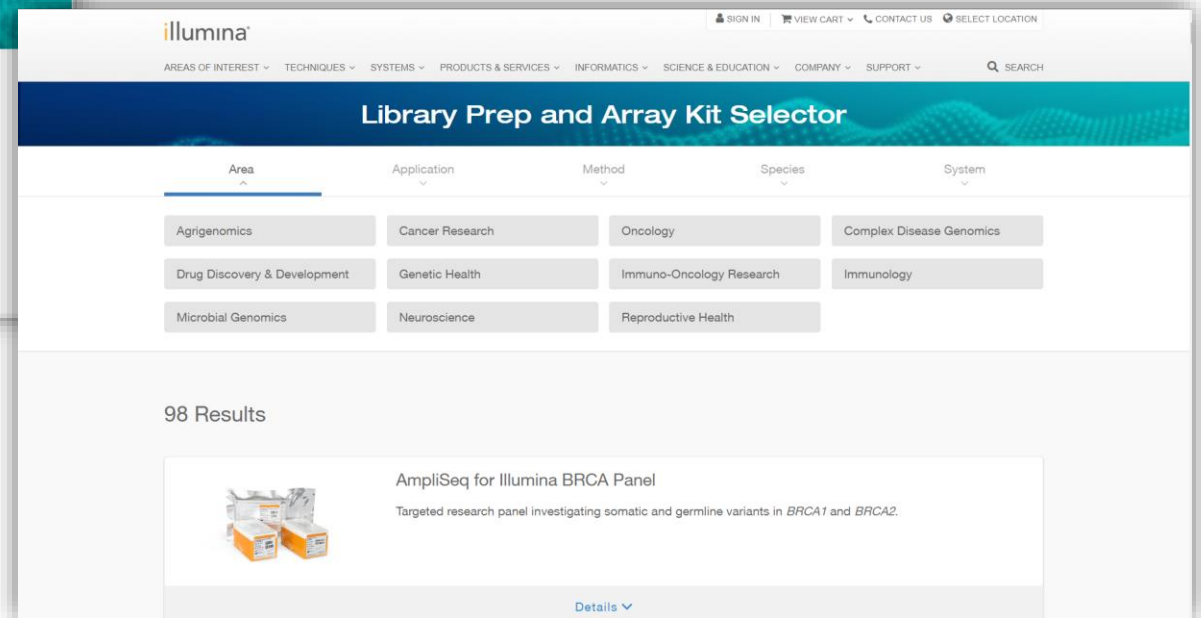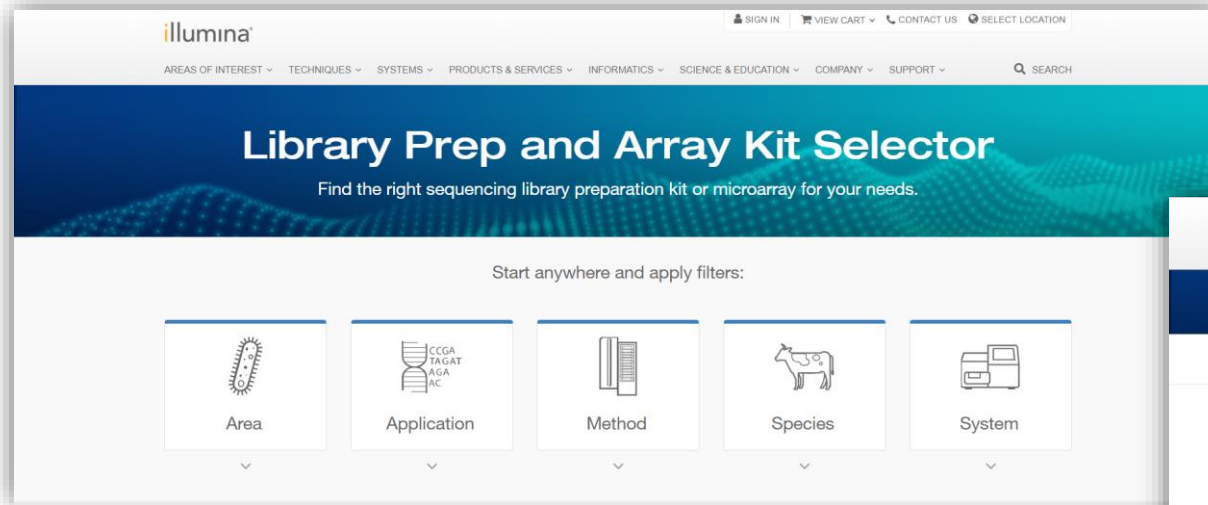

<http://www.illumina.com/library-prep-array-kit-selector.html>




MiSeq Sequencing Technology

Illumina Sequencing Workflow


1

Library Preparation

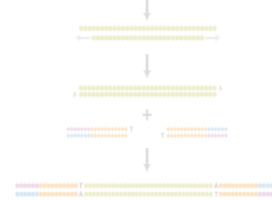


2

Cluster Generation

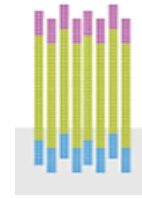

3

Sequencing

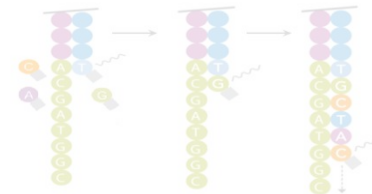


4

Data Analysis

```
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
```

Cluster Generation Overview

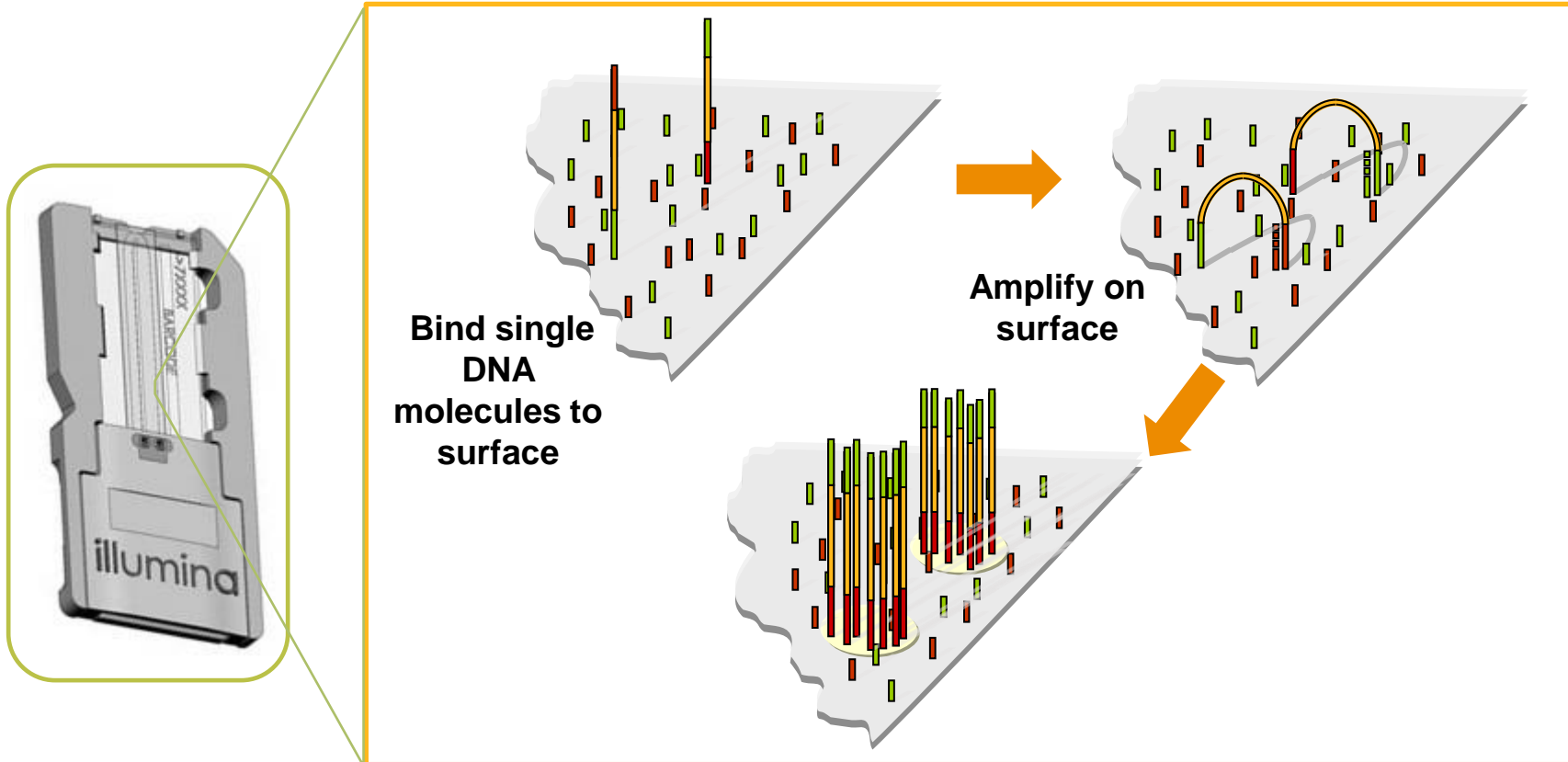

If the internet is available, click [here](#) to view a video:

Data Yield Increases as Number of Usable Cluster Increases


Cluster is unit of data production



Illumina Sequencing Workflow


1

Library Preparation

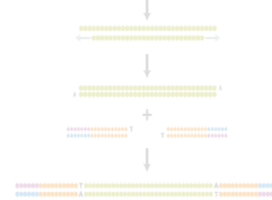


2

Cluster Generation

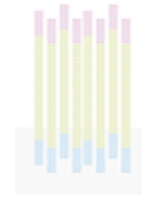

3

Sequencing

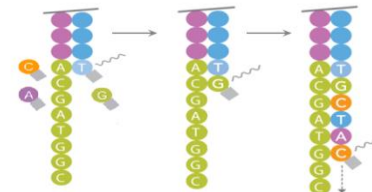


4

Data Analysis


```
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
```

Sequencing by Synthesis

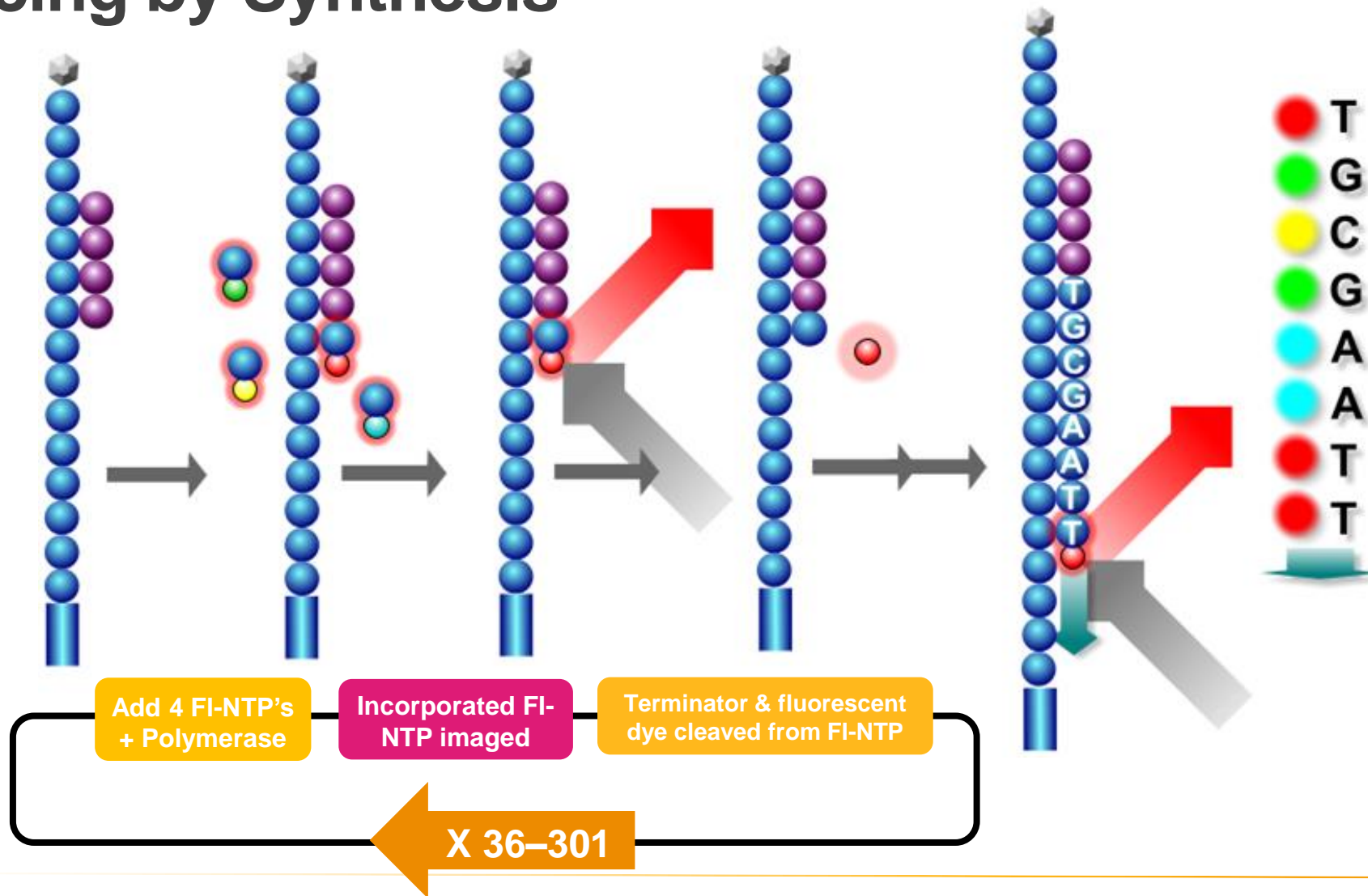

Imaging on the MiSeq

- Clusters are imaged using LEDs, cameras and filter combinations specific for each fluorescently labeled nucleotide

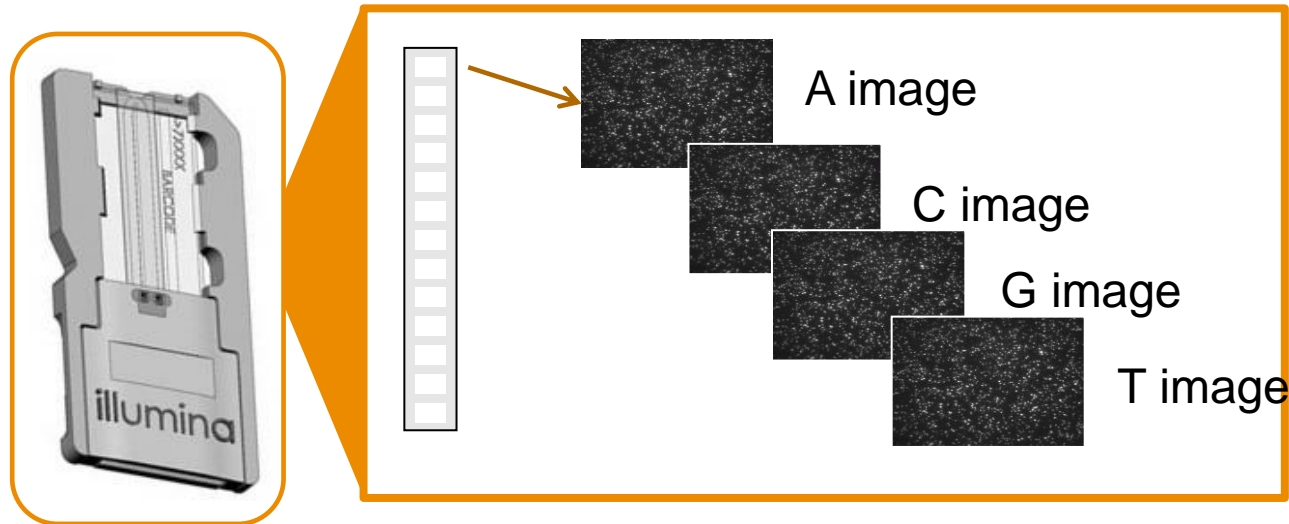


- After imaging is complete for one section (tile), the flow cell is moved to the next tile and the process is repeated

Illumina Sequencing Workflow


1

Library Preparation

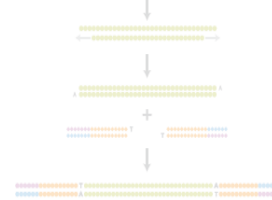


2

Cluster Generation

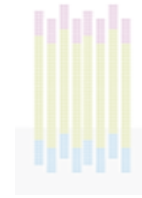

3

Sequencing

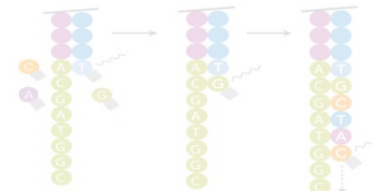

4

Data Analysis

```
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
STAAGGCTAGGTTTCATGCTA
```



SUBTITLE

Introduction to BaseSpace Sequence Hub



Analysis Tools

Diverse and growing set of tools allow for easier analysis of WGS data and new insights

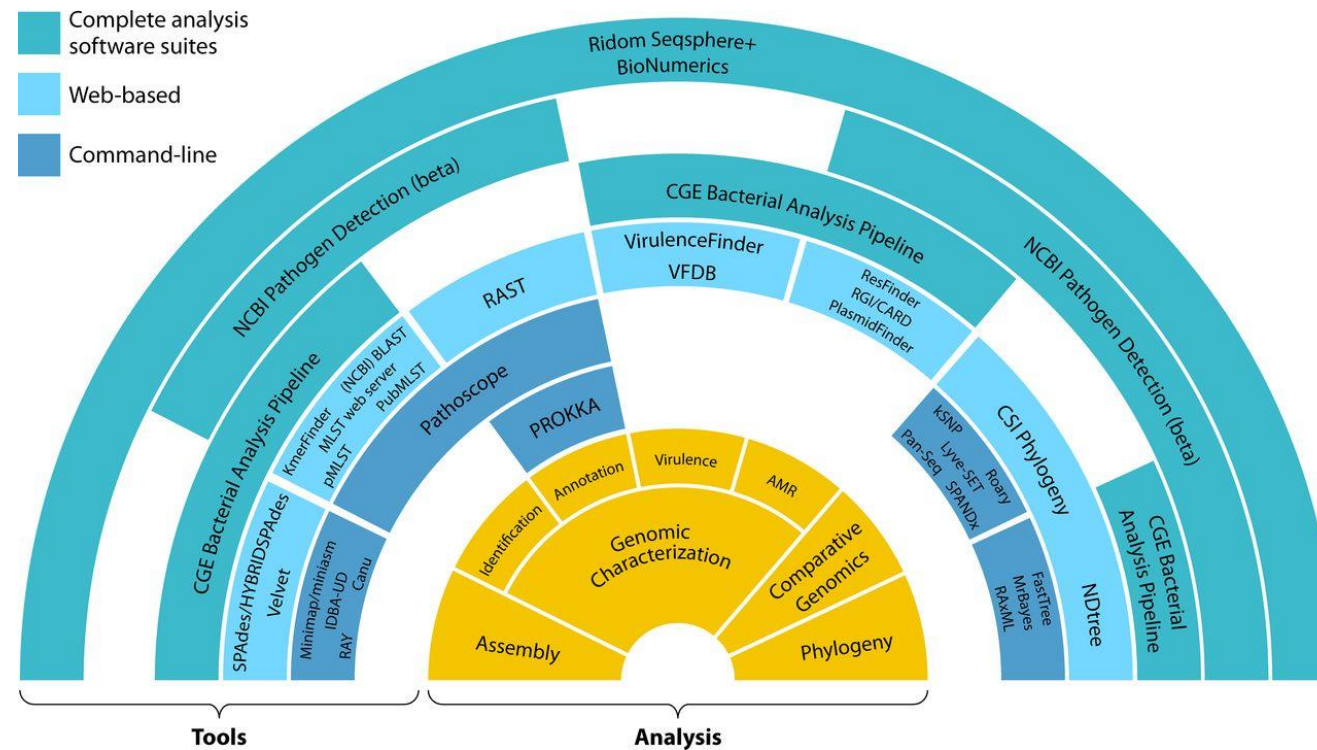

Clinical Microbiology Reviews Aug 2017, 30 (4) 1015-1063; DOI: 10.1128/CMR.00016-17



Vision

The destination for genomic investigators to **manage** and **analyze** genomic data, collaborate, and drive towards biological insights.

For Research Use Only. Not for use in diagnostic procedures.

BaseSpace Sequence Hub at-a-Glance*

7,500
Instruments
Connected

113K
Uploaded
Runs

95
Published
Apps

>20 PB
Data
Stored

40,000
Unique
Global
Users

**Instances in
USA, China,
EU, Aus,
Canada**

*As of 2019

For Research Use Only. Not for use in diagnostic procedures.

Accessing BaseSpace Sequence Hub

- US Instance
 - AWS Virginia
 - More compute power
 - <https://basespace.illumina.com>
- AU Instance
 - AWS Sydney
 - Data Stores within Australia
 - <https://aps2.sh.basespace.illumina.com>
- 2 instances share the same set of credentials
 - Data cannot be shared/transferred across regions

BaseSpace™ Sequence Hub

Enables anyone to analyze their NGS data




BaseSpace™ Sequence Hub

Genomics computing environment for NGS data

Monitor run quality,
share data


Upload directly
from instrument

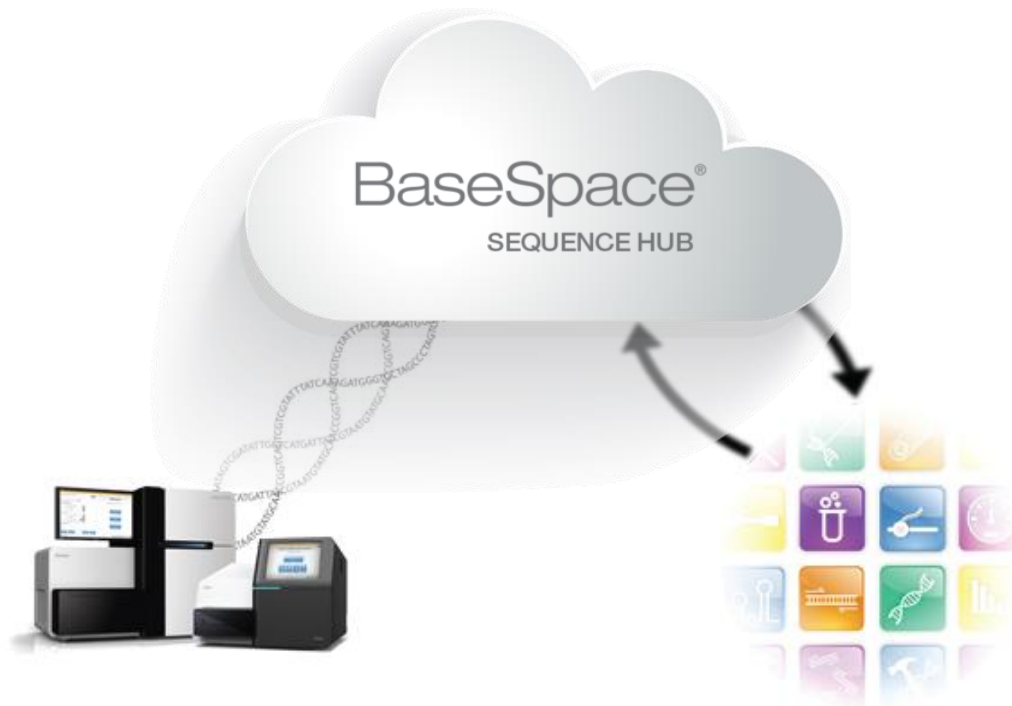

Share and
Collaborate




Broad Range of Sequence Hub Applications for Secondary Analysis

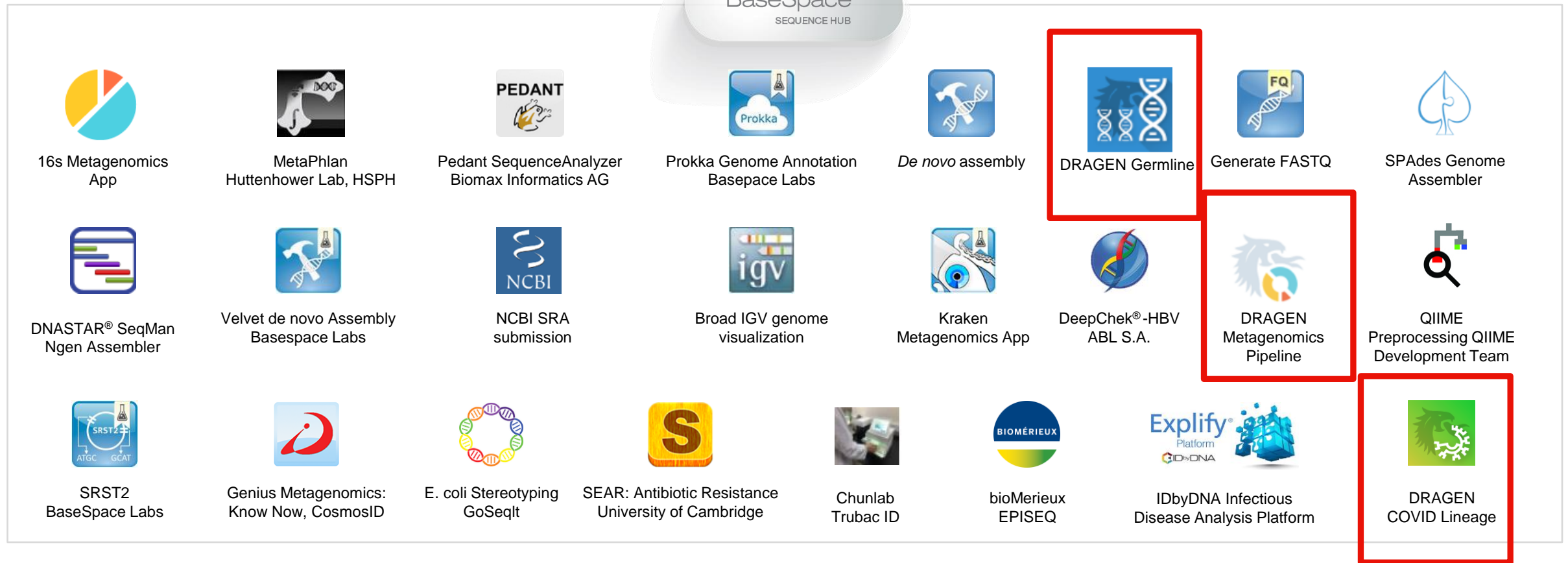

DRAGEN™ is Hardware-Accelerated Secondary Analysis


Dynamic Read Analysis for GENomics

DRAGEN Software Pipelines

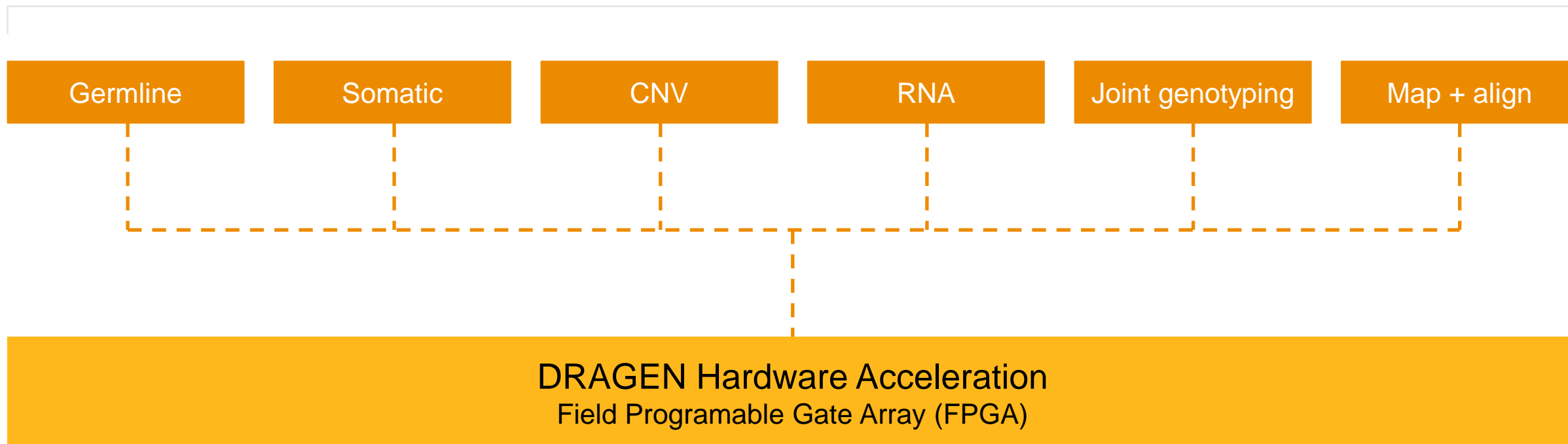

Flexible Data Analysis

On-premise

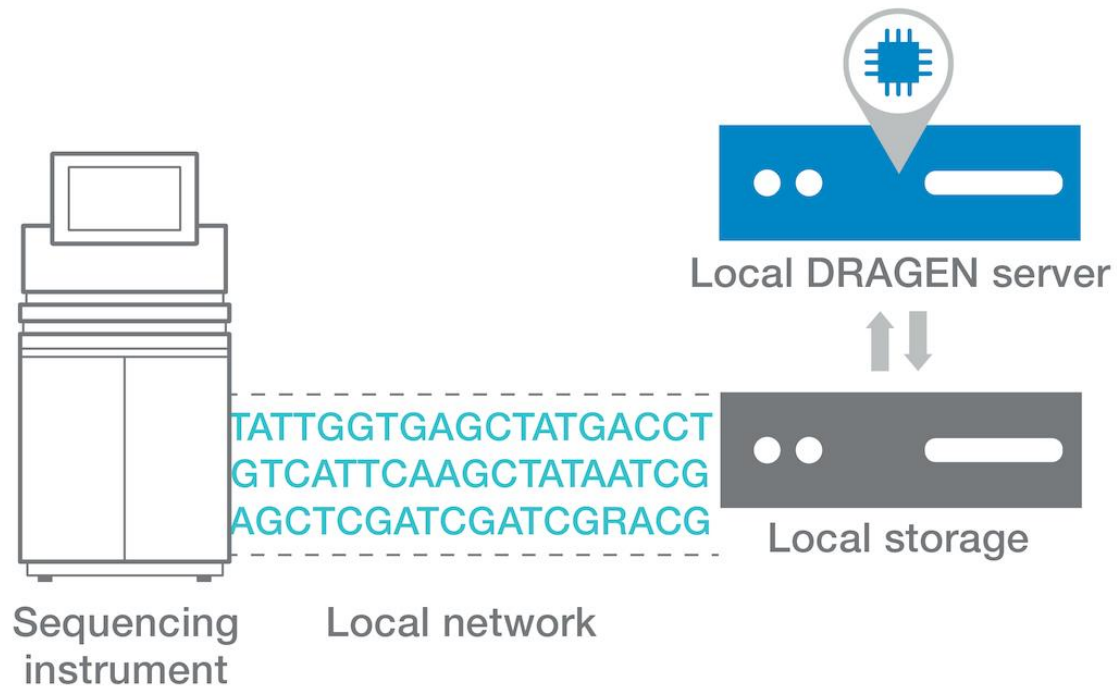

BaseSpace™ Sequence Hub



DRAGEN on BaseSpace Sequence Hub

Accurate, rapid secondary analysis in an easy-to-use, cloud-based environment



Accurate, Fast Analysis



Simple Workflow



Low-cost, Scalable Platform



Secure, Compliant Environment

Available Pipelines


- ✓ DRAGEN Germline Pipeline
- ✓ DRAGEN Somatic Pipeline
- ✓ DRAGEN Enrichment Pipeline
- ✓ DRAGEN RNA Pipeline
- ✓ DRAGEN Joint Genotyping Pipeline
- ✓ DRAGEN Methylation Pipeline
- ✓ DRAGEN Reference Builder

Bacterial Whole Genome Alignment and Variant Calling




Bacterial Whole genome resequencing


- Microbial whole-genome resequencing involves sequencing the entire genome of a bacteria, virus, or other microbe, and comparing the sequence to that of a known reference.
- The Germline (resequencing) workflow compares the DNA sequence in the samples against a reference genome and identifies any variants (SNPs or indels) relative to the reference sequence.

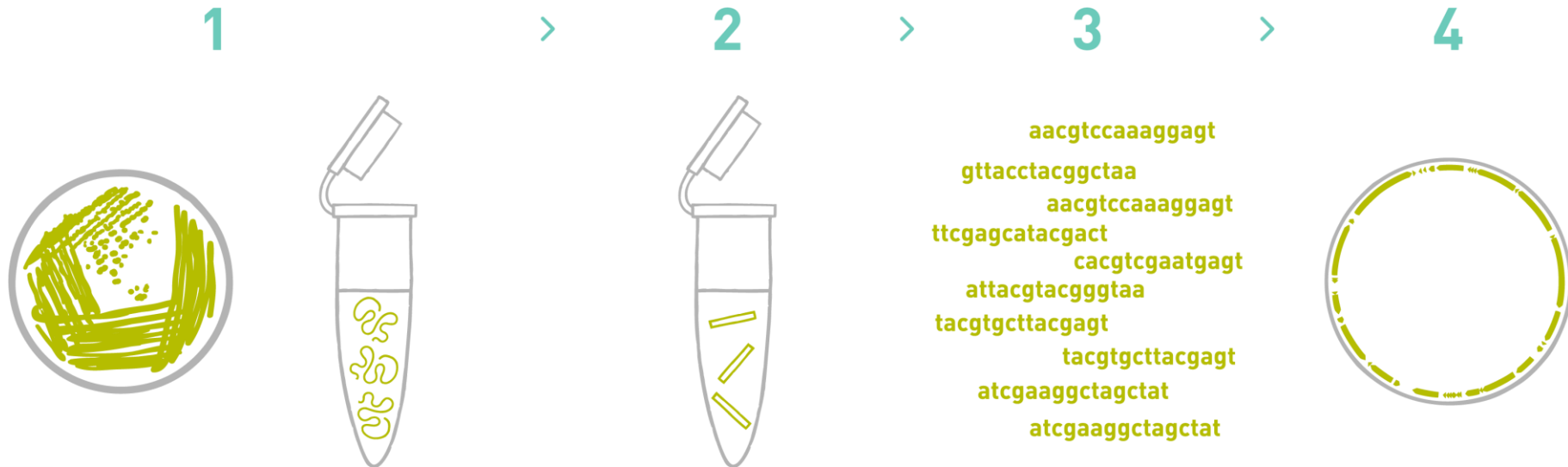


Whole Genome Resequencing

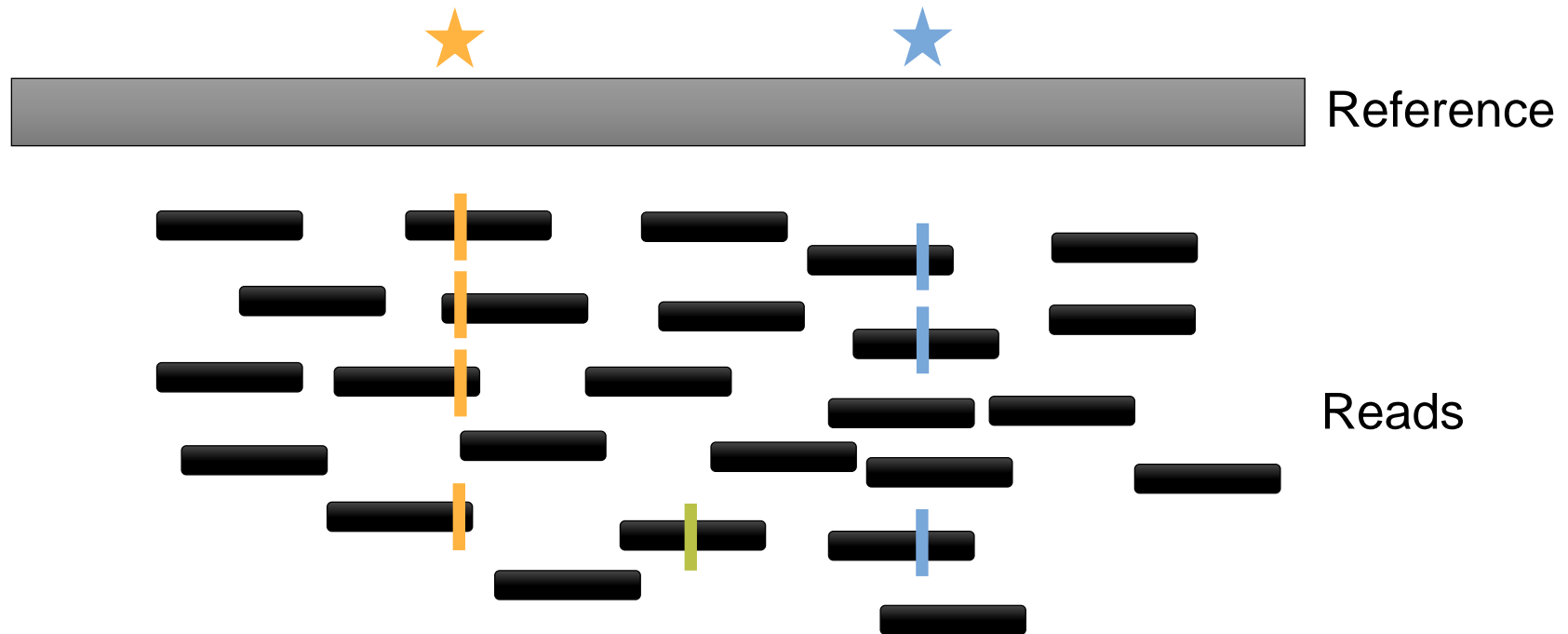



What is variant calling?

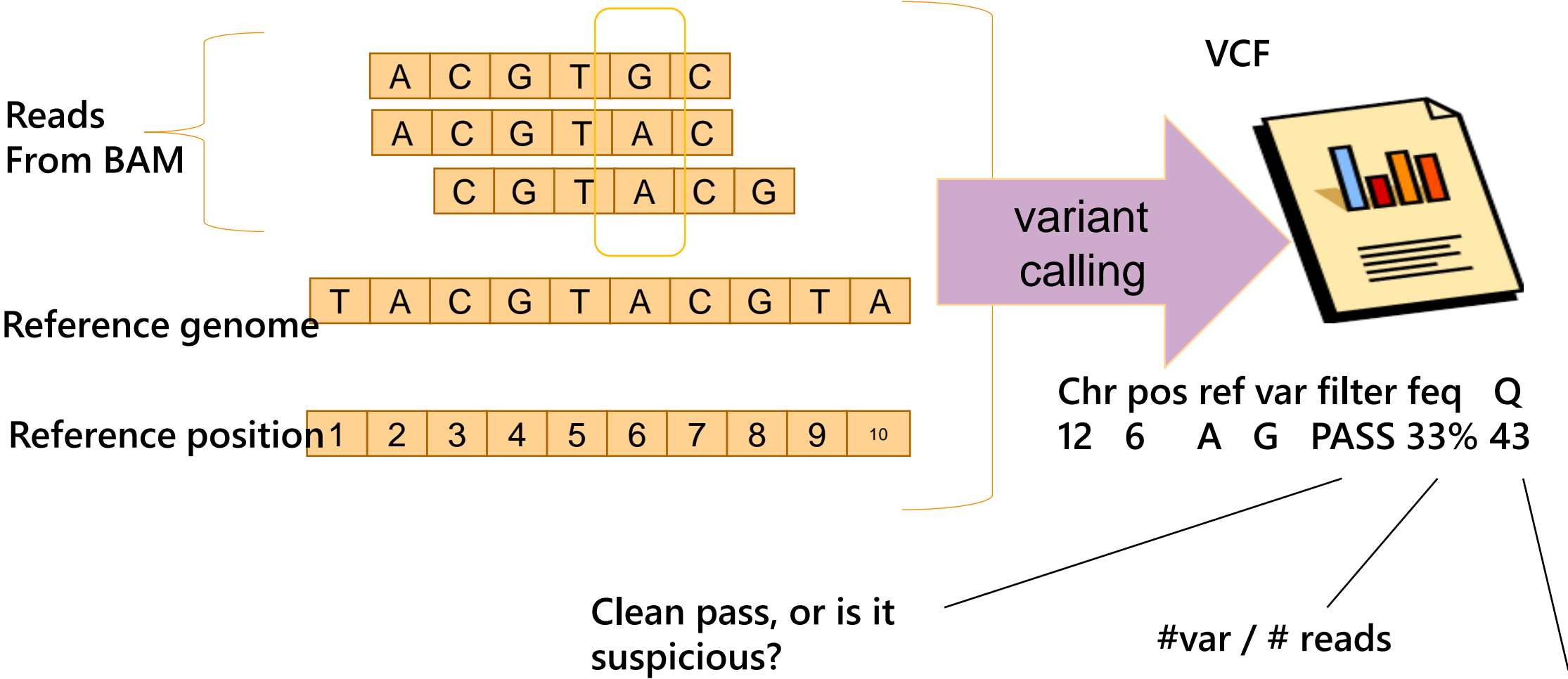


DRAGEN Germline Pipeline

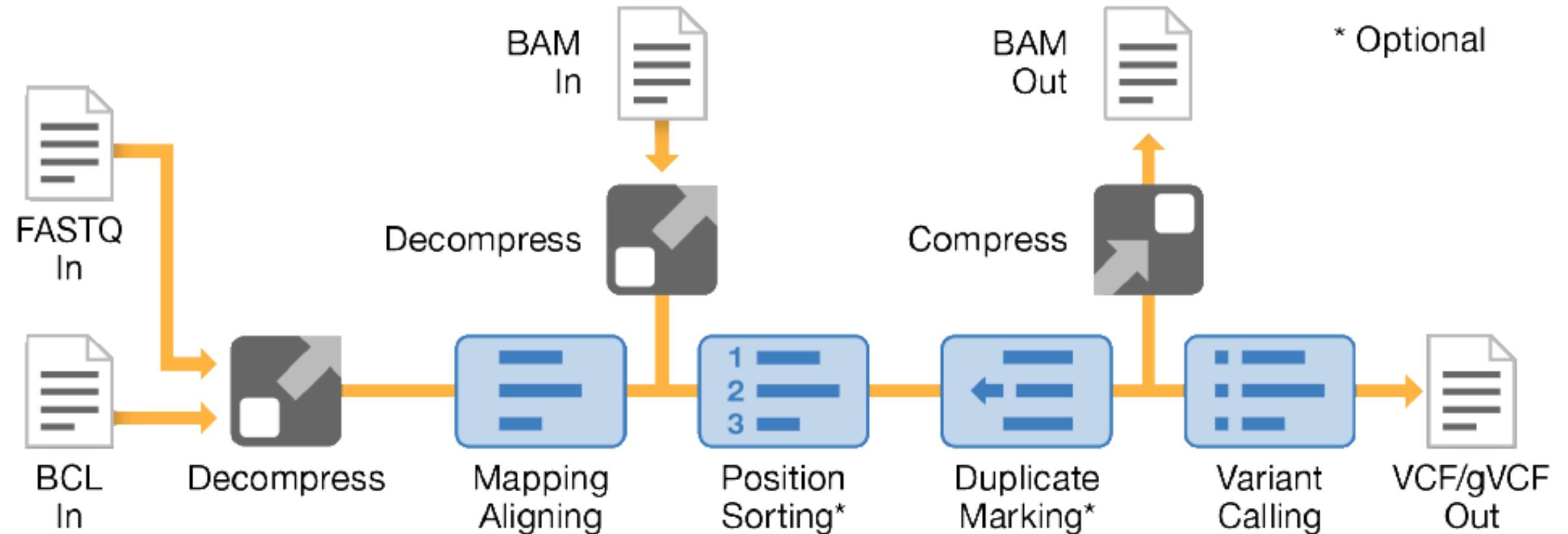


Mapping & Aligning—Germline Pipeline Example

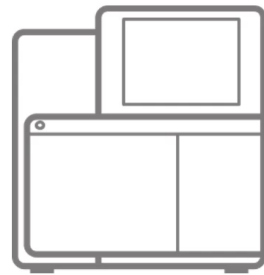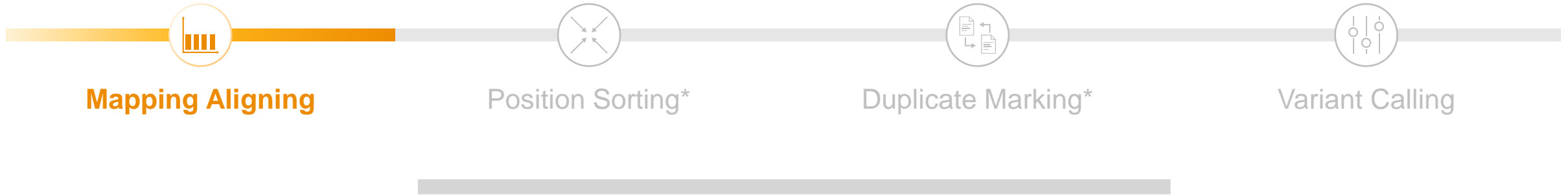


Position Sorting—Germline Pipeline Example

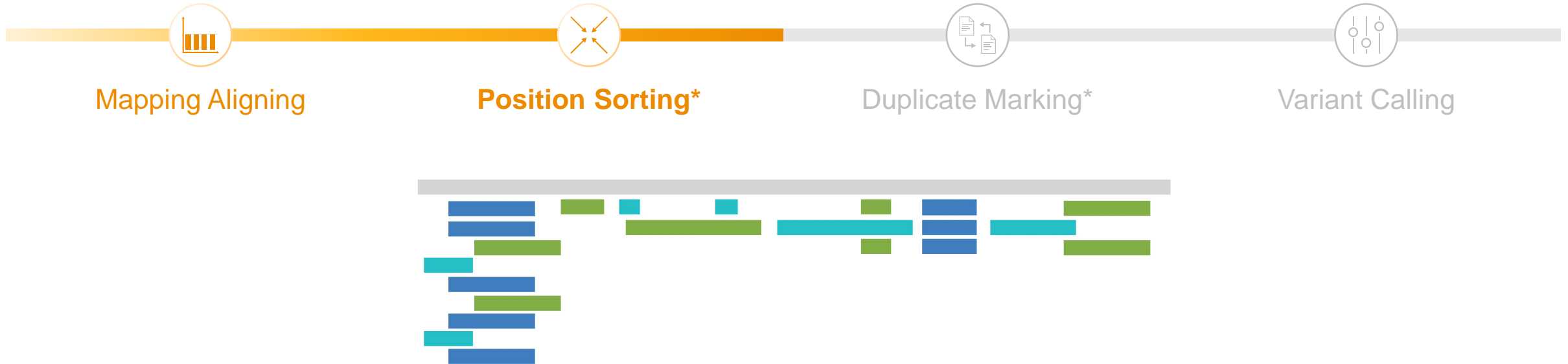

Duplicate Marking—Germline Pipeline Example

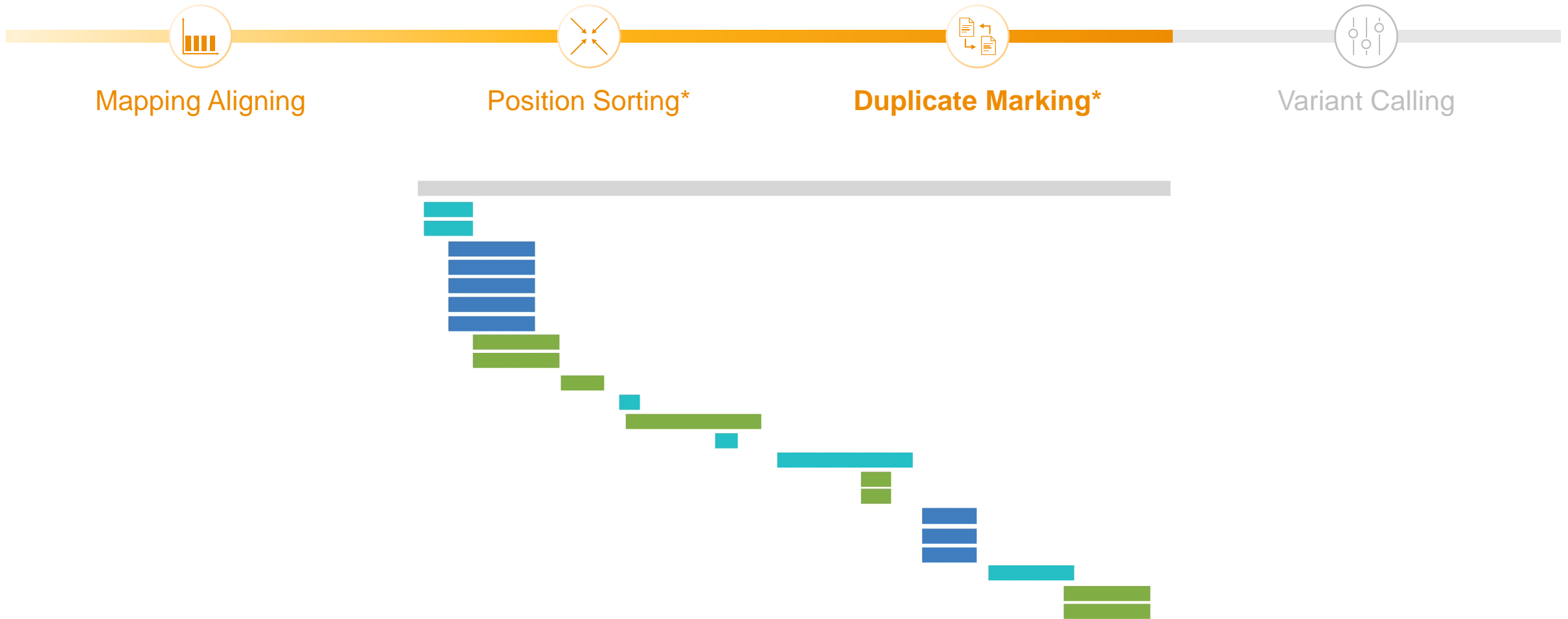

Variant Calling—Germline Pipeline Example

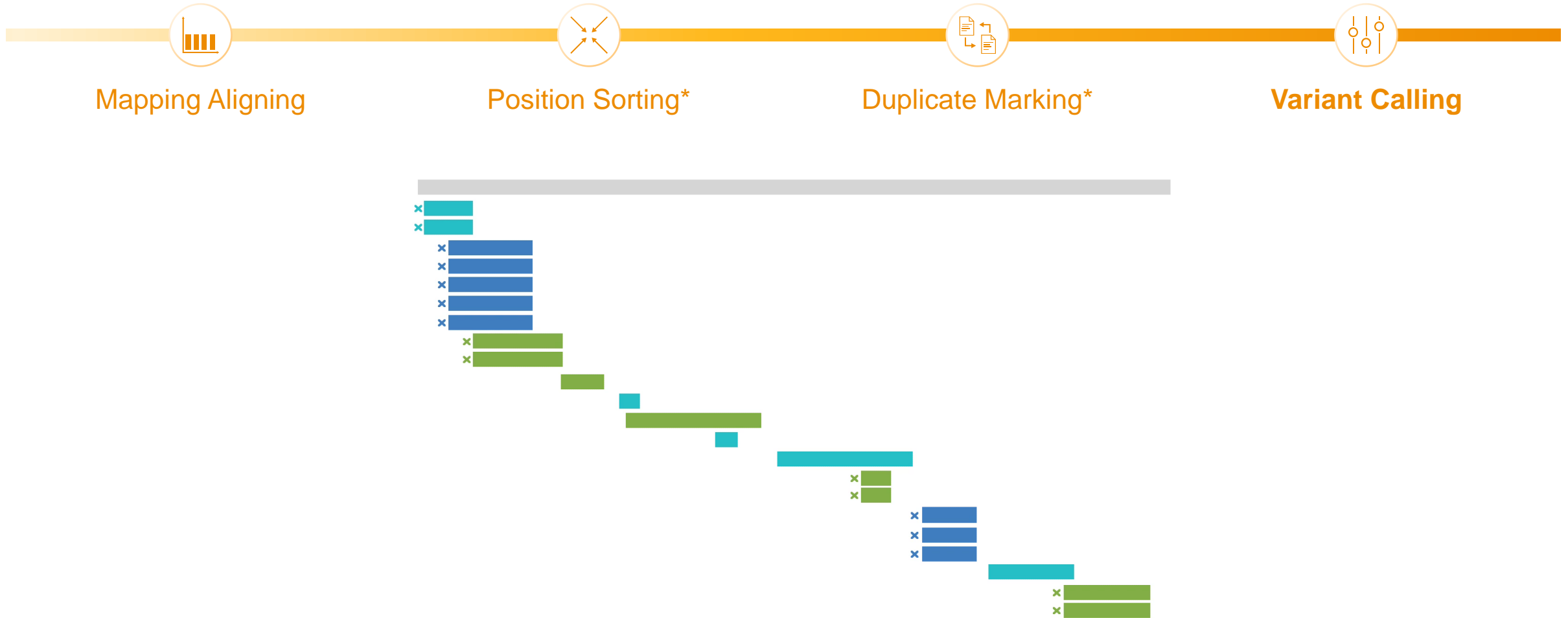

DRAGEN Automatic QC Metrics Reporting: Mapper

Removes the need to run downstream tools for QC

- ✓ Number of samples
- ✓ Reads Processed
- ✓ Total
- ✓ Biallelic
- ✓ Multiallelic
- ✓ SNPs
- ✓ INDELs
- ✓ MNPs
- ✓ SNP Transitions
- ✓ SNP Transversions
- ✓ Ti/Tv ratio

- ✓ Heterozygous
- ✓ Homozygous
- ✓ Het/Hom ratio
- ✓ In dbSNP
- ✓ Novel
- ✓ Total
- ✓ Biallelic
- ✓ Multiallelic
- ✓ SNPs
- ✓ INDELs
- ✓ MNPs

DRAGEN Automatic QC Metrics Reporting: Variant Caller

- Total input reads
- Number of duplicate reads (marked not removed)
- Number of unique reads
- Reads with mate sequenced
- Reads without mate sequenced
- QC-failed reads
- Mapped reads
- Number of unique & mapped reads (excl. dups)
- Unmapped reads
- Singleton reads (itself mapped; mate unmapped)
- Paired reads (itself & mate mapped)
- Properly paired reads
- Not properly paired reads (discordant)
- Reads with MAPQ [40:inf)
- Reads with MAPQ [30:40)
- Reads with MAPQ [20:30)
- Reads with MAPQ [10:20)
- Total reads in RG
- Supplementary (chimeric) alignments
- Average sequenced coverage over genome
- Total alignments
- Secondary alignments
- Supplementary (chimeric) alignments
- Estimated read length
- Bases in reference genome
- Bases in target bed [% of genome]
- Average sequenced coverage over genome
- Average alignment coverage over genome
- PCT of genome with coverage [40x:inf)
- PCT of genome with coverage [30x:40x)
- PCT of genome with coverage [20x:30x)
- PCT of genome with coverage [10x:20x)
- PCT of genome with coverage [5x:10x)
- PCT of genome with coverage [2x: 5x)
- PCT of genome with coverage [1x: 2x)
- PCT of genome with coverage [0x: 1x)
- DRAGEN mapping rate [mil. reads/second]
- Secondary alignments
- Estimated read length
- Insert length: mean
- Number of duplicate reads (marked)
- Number of unique reads
- Reads with mate sequenced
- Reads without mate sequenced
- QC-failed reads
- Mapped reads
- Number of unique & mapped reads (excl. dups)
- Unmapped reads
- Singleton reads (itself mapped; mate unmapped)
- Paired reads (itself & mate mapped)
- Properly paired reads
- Not properly paired reads (discordant)
- Reads with MAPQ [40:inf)
- Reads with MAPQ [30:40)
- Reads with MAPQ [20:30)
- Reads with MAPQ [10:20)
- Reads with MAPQ [0:10)
- Total alignments
- Secondary alignments
- Insert length: standard deviation

WGS Tools Help Genomic Epidemiology Characterize the Flow of Resistomes Between Animals and Humans



Bioinformatics tools available to detect Antimicrobial Resistance (AMR) determinants

- Public Health England's GeneFinder
- NCBI's AMRFinder
- KmerResistance
- SRST2
- ResFinder



WGS can predict AMR and the minimal inhibiting concentration (MIC) of an antimicrobial, applying machine or deep learning to genome sequence data^[1].



FDA, the Centers for Disease Control and Prevention (CDC), and the U.S. Department of Agriculture (USDA) established the National Antimicrobial Resistance Monitoring System (NARMS) to systematically sequence pathogens from food-animals for AMR surveillance^[2].

1. Scientific Report. (2018) 8:421–18972. doi: 10.1038/s41598-017-18972-w
2. Foodborne Pathog Dis. (2017) 14:545–57. doi: 10.1089/fpd.2017.2283


Culture-based Whole-Genome Sequencing (WGS) in Action

Zoonotic transfer of *Clostridium difficile* harboring antimicrobial resistance between farm animals and humans

“The phylogeny demonstrates clear mixing of European and North American strains, indicating multiple transmission events between continents, and mixing of human and animal strains, indicating multiple transmissions events between these hosts.”

Sequenced on Illumina HiSeq™ Sequencer

C. W. Knetsch et al. Journal of Clinical Microbiology Feb 2018, 56 (3) e01384-17; DOI: 10.1128/JCM.01384-17

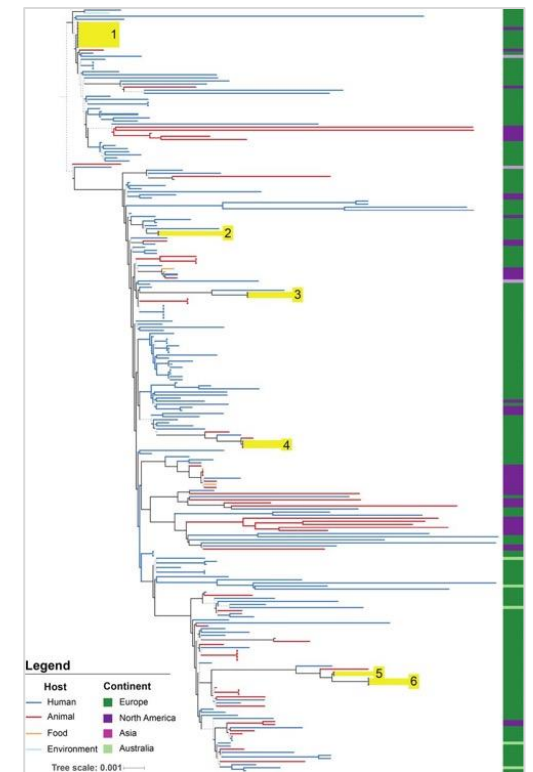

|| Microbial Genomics on the MiSeq

Advantage of Using Illumina Technology in Microbial Studies

- Sequencing the 5 Mb genome of *Escherichia coli* with Illumina technology can be done in one day at a fraction of the cost of traditional capillary electrophoresis methods
- Microbial genomes are small and data analysis is relatively simple
- NGS has the ability to measure the changes in the genome with prior knowledge
- Single-base resolution allows tracking of microbial adaptation over a short period of time



Microbial Genomics Resources

illumina

Search Illumina.com

Products Learn Company Support Recommended Links

MICROBIAL GENOMICS Overview Sequencing Methods Microbiome Analysis Infectious Diseases Environmental Metagenomics More

QUESTIONS

Accelerate microbial breakthroughs

Areas of Interest:
Microbial Genomics

Illumina next-generation sequencing takes you inside microbiology

Putting the big picture into focus by revealing the smallest of details

Characterize un-culturable organisms. Discover entirely new viruses. Develop new strategies to control outbreaks. Monitor host-pathogen interactions.


Next-generation sequencing (NGS) is opening new doors in microbial genomics, revealing fresh insight into how microbes impact humans and the environment.

- See how NGS enables a broad range of microbial studies. [View Video](#).
- Explore how Dr. Jennifer Gardy is using NGS to provide new clues to help prevent and control disease outbreaks. [Read Article \(PDF\)](#).
- Learn how Dr. Ramunas Stepanauskas uses single-cell sequencing on Illumina systems to uncover hidden diversity and biogeographic variability in marine

Microbiology Methods Guide

All the information you need, from BeadChips to library preparation to sequencer selection and analysis. Select the best tools for your lab.

[Access Guide](#)

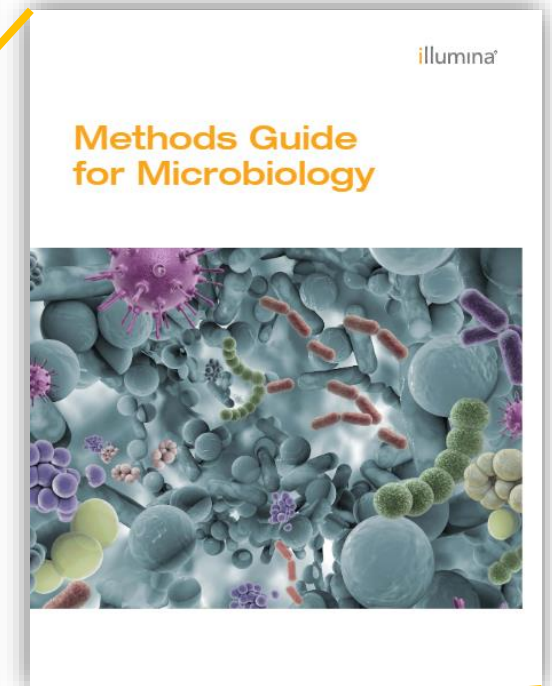

<https://www.illumina.com/areas-of-interest/microbiology.html>

|| Thank you

For any additional technical queries please
contact Illumina Tech Support

techsupport@illumina.com or 1800.775.688